

Cells and Tissue DNA Isolation Kit (Magnetic Bead System) – 50 Preps Product # 59100

Product Insert

Norgen's Cells and Tissue DNA Isolation Kit (Magnetic Bead System) provides a fast and reproducible method for isolating genomic DNA from cultured cells as well as various tissue samples, bodily fluids and nasal or throat swabs. The purified DNA is of the highest quality, and is compatible with a number of downstream research applications including PCR, NGS and microarray analysis. Using this kit, DNA can be extracted either manually or by using automation platform IsoPureTM. The protocol can also be easily modified by making minor changes to be able to run on other Magnetic Bead-based automation platforms such as KingFisherTM Flex 96 and Hamilton MagEx STAR platforms.

Norgen's Purification Technology

Purification is based on the use of magnetic beads that bind DNA under optimized binding conditions. Norgen's Cells and Tissue DNA Isolation Kit (Magnetic Bead System) allows for the isolation of genomic DNA from various types of animal tissues or cell samples. Different sample lysate preparation steps are applied, including a Proteinase K treatment. Magnetic Beads A and ethanol are added to the clean supernatant, and the resulting solution is placed on the magnetic beads, while most of the proteins will be removed in the supernatant. The bound DNA is then washed with Solution WN and 70% ethanol in order to remove any remaining impurities and the purified total DNA is eluted with the Elution Buffer B. The purified DNA can be used in a number of downstream applications.

Specifications

Kit Specifications			
Number of Preps	50		
Maximum Cells and Tissue Input	20 mg of animal tissue 3 x 10 ⁶ cells		
Average Yield *	8-10 μg (20 mg of animal tissue) 8-12 μg (3 x 10 ⁶ cells)		
Average Purity (OD260/280)	1.8 – 1.9		
Time to Complete 50 Purifications (automated)	15 minutes (hands-on time)		
Time to Complete 50 Purifications (manual)40 minutes (hands-on time)			

* Average DNA yield will vary depending on the type of samples.

Kit Components

Component	Product #59100 (50 samples)
Lysis Buffer B	20 mL
Proteinase K in Storage Buffer	1.2 mL
Magnetic Beads A	2.2 mL
Solution WN	18 mL
Elution Buffer B	15 mL
Elution Tubes (1.7 mL)	50
Product Insert	1

Storage Conditions and Product Stability

All solutions should be kept tightly sealed and stored at room temperature. This kit is stable for 1 year after the date of shipment. The kit contains a ready-to-use Proteinase K, which is dissolved in a specially prepared storage buffer. The buffered Proteinase K is stable for up to 1 year after the date of shipment when stored at room temperature.

Advantages

- Isolate genomic DNA from cultured cells as well as various tissue types
- Fast, reproducible and easy processing using a magnetic bead system
- Isolate high quality genomic DNA
- Recovered genomic DNA is compatible with various downstream applications
- Compatible with automation robotic system

Precautions and Disclaimers

This kit is designed for research purposes only. It is not intended for human or diagnostic use. Ensure that a suitable lab coat, disposable gloves and protective goggles are worn when working with chemicals. For more information, please consult the appropriate Material Safety Data Sheets (MSDSs). These are available as convenient PDF files online at *www.norgenbiotek.com*.

Solution WN contains guanidinium salts and should be handled with care. Guanidinium salts forms highly reactive compounds when combined with bleach, thus care must be taken to properly dispose of any of these solutions.

Cells and Tissue of all human and animal subjects is considered potentially infectious. All necessary precautions recommended by the appropriate authorities in the country of use should be taken when working with these samples.

Customer-Supplied Reagents and Equipment

- Magnetic bead separation rack
- IsoPure[™] 96 or IsoPure[™] Mini Purification System (For automated Isolation)
- Micropipettors
- Microcentrifuge tube
- 70% ethanol (prepare fresh)
- 96 100% ethanol
- Temperature adjustable (37°C, 55°C, 65°C) incubator(s)
- Nuclease-free water
- Phosphate Buffered Saline (PBS) (for Cultured Animal Cells protocol)
- Cell Disruption Tools such as mortar and pestle, rotor-stator homogenizer (for Animal Tissue protocol)
- Syringe with a 22G needle (for Animal Tissue protocol)
- Orbital Shaker (Plates) or Tube Revolver Rotator (Tubes)
- RNase A (optional)


Procedure

Notes prior to use:

- The procedure provides information and steps to follow to extract DNA using automated as well as manual method.
- In each method, the steps for preparing the lysate are different depending on the starting material (Section 1). However, the subsequent steps are the same in all cases (Section 2).
- Please ensure that the correct procedure for preparing the lysate from your starting material is followed.
- Ensure that all solutions are at room temperature prior to use, and that no precipitates have formed. If necessary, warm the solutions and mix well until the solutions become clear again.
- Always vortex the Magnetic Beads A before use.
- Always vortex the Proteinase K before use.
- Preheat the incubator(s) according to the temperatures required (37°C or 55°C or 65°C). For automation procedures, the automation platform that is being used may have a heat block to provide appropriate temperatures for incubation and thus an incubator might not be required.
- Prepare a working concentration of the Solution WN by adding 24 mL of 96 100 % ethanol (provided by the user) to the supplied bottle containing the concentrated Solution WN. This will give a final volume of 42 mL. The label on the bottle has a box that may be checked to indicate that the ethanol has been added.

Flow Chart

Procedure for using Norgen's Cells and Tissue DNA Isolation Kit (Magnetic Bead System) – IsoPure[™] 96

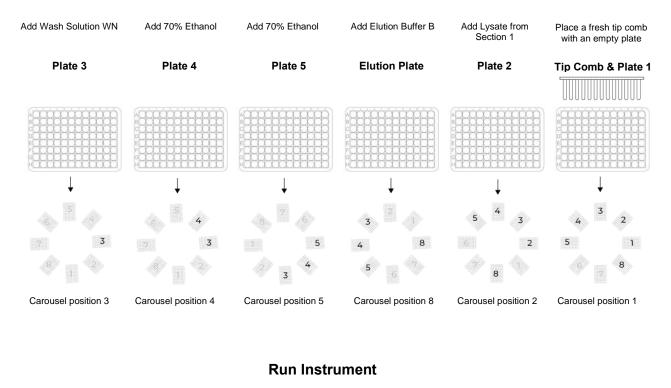




Add ethanol & bead mix to Plate 2 when prompted Press "OK" after placing plate 2 back at position 2

Residence.

Run Complete. Take Elution plate and store DNA.

Flow Chart

Procedure for using Norgen's Cells and Tissue DNA Isolation Kit (Magnetic Bead System) – IsoPure[™] Mini

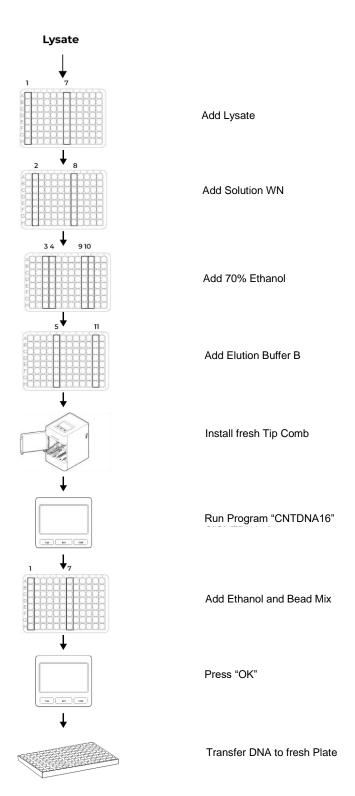

I. Genomic DNA Isolation using Automation Platforms

Procedure for IsoPure[™] Mini or IsoPure[™] 96 Automation System

The procedure outlined below is for extraction of DNA using automation platforms IsoPure[™] 96 (Cat. AP1096) and IsoPure[™] Mini (Cat. AP1016). This procedure can easily be adapted to other automation platforms such as KingFisher[™] Flex 96 and Hamilton MagEx STAR platforms.

Set-up the IsoPure[™] Plates

Notes Prior to Use

- Only use plates that are compatible with the instrument that is being used. Use 96 deep well plates for IsoPure 96 (Cat. AP1096-DWP), 96 well elution plates for IsoPure[™] (Cat. AP1096-ELP) and Tip Comb for IsoPure 96 (Cat. AP1096-TC) while working with IsoPure[™] 96 Instruments. Use 96 deep well plates for IsoPure 96 (Cat. AP1096-DWP) and 8-place magnetic tip comb for IsoPure Mini (Cat. AP1016-MC) while working with IsoPure[™] Mini Instruments.
- The numbering of plates is based on where it should be placed on the carousel in the instrument. Use the carousel rotation keys on the instrument or select the carousel position on the screen to bring it to the right position as mentioned in the flow chart.
- IsoPure[™] 96 uses multiple plates for DNA extraction while IsoPure[™] Mini uses single plate for DNA extraction.
- For IsoPure[™] 96 most of the lysate preparation is performed in plate 2 while other plates are used for performing washing and elution steps.
- For IsoPure[™] Mini the lysate preparation is performed in column 1 and/or 7 while washing and elutions steps could be performed in other columns of the same plate.
- For IsoPure[™] 96, DNA elution step is always performed in plate 8.
- Prepare the plates and label them based on the sample type as mentioned in the tables below.

Table 1. Set-up for Animal Tissue DNA Extraction

Prepare plates as per the volumes mentioned in the table below. Detailed description of step-wise process can be found in **Section 1 and 2.**

Component	IsoPure™ Mini	IsoPure [™] 96	Location on Carousel	Amount
Tip Comb	In the instrument	Plate 1	1	Empty Plate
Lysate	Column 1 / Column 7	Plate 2	2	Variable
Wash WN	Column 2 / Column 8	Plate 3	3	500 µL
70 % Ethanol	Column 3 / Column 9 and Column 4 / Column 10	Plate 4 & Plate 5	4 & 5	500 µL
Elution Buffer B	Column 5 / Column 11	Elution Plate	8	50-100 μL

Table 2. Set-up for Cells and Viral Suspension DNA Extraction

Prepare plates as per the volumes mentioned in the table below. Detailed description of step-wise process can be found in **Section 1 and 2.**

Component	IsoPure [™] Mini	IsoPure [™] 96	Location on Carousel	Amount
Tip Comb	In the instrument	Plate 1	1	Empty Plate
Sample	Column 1 / Column 7	Plate 2	2	Variable
Proteinase K	Column 1 / Column 7	Plate 2	2	20 µL
Lysis Buffer B	Column 1 / Column 7	Plate 2	2	300 µL
Wash WN	Column 2 / Column 8	Plate 3	3	500 µL
70 % Ethanol	Column 3 / Column 9 and	Plate 4 & 5	4 & 5	500 µL
	Column 4 / Column 10			
Elution Buffer B	Column 5 / Column 11	Elution Plate	8	50-100 µL

Section 1. Lysate Preparation from Various Sample Types

1A. Lysate Preparation from Animal Tissues

Notes Prior to Use

- Fresh or frozen tissues may be used for the procedure.
- Tissues should be immediately frozen and stored at -20°C or -70°C. Tissues may be stored at -70°C for several months.
- It is recommended that no more than 20 mg of tissue be used.

1A. Lysate Preparation from Animal Tissues

- a) Excise up to 20 mg of tissue sample. Either frozen or fresh tissue may be used.
- **b)** Homogenize the sample into a fine powder in liquid nitrogen using a mortar and pestle, then transfer into a nuclease-free microfuge tube.
- c) Add 300 µL Lysis Buffer B to the tissue sample.

Note: The sample can also be homogenized using tools such as a rotor-stator homogenizer or a microfugesize pestle in the provided Lysis Solution.

- d) Use a syringe with a 22G needle to further homogenize the lysate by passing through the needle 10 times.
- e) Add 20 µL of Proteinase K (vortex before use) to the lysate. Mix well by gentle vortexing.
- f) Transfer the entire lysate to Plate 2 (IsoPure[™] 96) or Column 1/Column 7 (IsoPure[™] Mini).
- g) Proceed to Section 2: DNA Isolation from All Types of Lysates.

1B. Lysate Preparation from Cultured Animal Cells

Notes Prior to Use

- Cells grown in suspension or monolayer may be used.
- The maximum recommended input of cells is 3 x 10⁶. A hemocytometer can be used in conjunction with a microscope to count the number of cells. As a general guideline, a confluent 3.5 cm plate of HeLa cells will contain 10⁶ cells.
- Cell pellets can be stored at -70°C for later use or used directly in the procedure. Determine the number of cells present before freezing.

Frozen cell pellets should not be subjected to repeated cycles of freeze and thaw prior to beginning the protocol.

1B (i). Cell Lysate Preparation from Cells Growing in a Monolayer

- Add 300 µL Lysis Buffer B and 20 µL Proteinase K to all the wells in Plate 2 (IsoPure[™] 96) or Column 1 and/or Column 7 (IsoPure[™] Mini) where the sample will be added. For ease of use, prepare a mixture of Lysis Buffer B and Proteinase K. (For example, if DNA is to be extracted from 10 samples, mix 3 mL of Lysis Buffer B and 200µL of Proteinase K). Make 10% extra volume to avoid volume loss due to pipetting.
- b) Detach cells by standard trypsinization method or cell scraper. Transfer an appropriate amount to a 1.5 mL microfuge tube (not provided). The maximum recommended input of cells is 3 x 10⁶.
- c) Collect cells by centrifugation at no more than 200 x g (~2,000 RPM) for 10 minutes. Discard the supernatant.
- d) Add 100 µL of phosphate buffered saline (PBS) to the cell pellet. Mix and transfer the cells dissolved in PBS to Plate 2 (IsoPure[™] 96) and Column 1 and/or Column 7 (IsoPure[™] Mini).
- e) Proceed to Section 2: DNA Isolation from All Types of Lysates.

1B (ii). Cell Lysate Preparation from Cells Growing in Suspension and Lifted Cells

a) Add 300 µL Lysis Buffer B and 20 µL Proteinase K to all the wells in Plate 2 (IsoPure[™] 96) and Column 1 and/or Column 7 (IsoPure[™] Mini) where the sample will be added. For ease of use, prepare a mixture of

Lysis Buffer B and Proteinase K. (for example, if DNA is to be extracted from 10 samples, mix 3 mL of **Lysis Buffer B** and 200 µL of **Proteinase K**). Make 10% extra volume to avoid volume loss due to pipetting.

- **b)** Transfer an appropriate amount of cells to a 1.5 mL microfuge tube (not provided). The maximum recommended input of cells is 3 x 10⁶.
- c) Collect cells by centrifugation at no more than 200 x g (~2,000 RPM) for 10 minutes. Discard the supernatant.
- d) Add 100 µL of phosphate buffered saline (PBS) to the cell pellet. Mix and transfer the cells dissolved in PBS to Plate 2 (IsoPure[™] 96) or Column 1 and/or Column 7 (IsoPure[™] Mini).
- e) Proceed to Section 2: DNA Isolation from All Types of Lysates.

1C. Lysate Preparation for Viral DNA

Notes Prior to Use

- For the isolation of integrated viral DNA, follow **Section 1A** if the starting material is animal tissue and follow **Section 1B** if the starting material is cell culture.
- For the isolation of DNA from free viral particles, follow the protocol provided below.
- Up to 150 µL of viral suspension can be processed.
- Fresh samples are recommended. Frozen samples may be used, however the yield of genomic DNA may be decreased.

1C. Lysate Preparation for Viral DNA

- a) Add 300 µL Lysis Buffer B and 20 µL Proteinase K to all the wells in Plate 2 (IsoPure[™] 96) and Column 1 and/or Column 7 (IsoPure[™] Mini) where the sample will be added. For ease of use, prepare a mixture of Lysis Buffer B and Proteinase K. (for example, if DNA is to be extracted from 10 samples, mix 3 mL of Lysis Buffer B and 200 µL of Proteinase K). Make 10% extra volume to avoid volume loss due to pipetting.
- b) Transfer up to 150 µL of viral suspension to Plate 2 (IsoPure[™] 96) and Column 1 and/or Column 7 (IsoPure[™] Mini).
- c) Proceed to Section 2: DNA Isolation from All Types of Lysates.

Section 2: DNA Isolation from All Types of Lysates

Table 3. 96-100% ethanol and magnetic bead volume to be added to lysate of each sample type.

Sample Type	Ethanol Volume (µL)	Magnetic Bead Volume (µL)	Total Volume (μL)
Tissue	300	40	340
Cells (Monolayer)	400	40	440
Cells (Suspension)	400	40	440
Viral Suspension	450	40	490

2.1 DNA Isolation using IsoPure[™] 96 System

- a) Setup all the plates according to **Table 1 or 2** as per the sample being isolated, add 500 µL of **Solution WN** into **Plate 3** and keep it at **position 3**.
- b) Add 500 µL freshly prepared 70% ethanol to plate 4 & 5 and keep it at position 4 & 5 respectively.
- c) Add 50-100 µL of Elution Buffer B to the elution plate and place it at position 8.
- d) Place sample plate (Plate 2) in the machine at position 2.
- e) Place a clean tip comb in an empty plate (Plate 1) and keep it at position 1.
- f) After setting up, run the program CNTDNA96.
- g) Prepare a mixture of 96-100% ethanol and Magnetic Beads A as per Table 3. For ease of pipetting, it is recommended to prepare a bulk mix containing 96-100% ethanol and Magnetic Beads A. We recommend a 10% overage to compensate for pipetting errors. Make sure to vortex the bulk mix before adding it to the lysate. It is recommended to prepare the bulk fresh every time a new DNA isolation is performed.

h) After 1 hour of incubation, the instrument will pause and will allow the user to add mixture of ethanol and magnetic beads prepared in step g. Open the cover and remove the plate from the instrument and add the prepared mixture from step g to Plate 2.

Optional RNase A Treatment:

If RNA-free genomic DNA is required, add the equivalent of 10 KUnitz of RNase A (not to exceed 20 μ L) to the lysate after addition of **96-100% ethanol** and **Magnetic Beads A** mixture to **Plate 2** during step g. If RNase treatment is performed, make sure to use **CNTDNA96R** (IsoPureTM 96) program instead of **CNTDNA96** to include an incubation step of 37°C for 15 minutes.

- i) Put the plate back into the instrument and press OK. The instrument will continue to perform the extraction.
- j) The instrument will prompt after the procedure is complete. Remove the elution plate from the instrument and the DNA is now ready for further downstream processing.
- k) The purified DNA sample may be stored at 4°C for a few days. It is recommended that samples be placed at -20°C for long-term storage.

2.2 DNA Isolation using IsoPure[™] Mini System

- a. Setup the plates according to **Table 1 or 2** as per the sample being isolated, add 500 µL of **Solution WN** into **Column 2/8**.
- b. Add 500 μL freshly prepared 70% ethanol to Column 3 and/or Column 9 and Column 4 and/or Column 10.
- c. Add 50-100 µL of Elution Buffer B to Column 5 and/or Column 11.
- d. Place a clean tip comb into the IsoPure[™] Mini.
- e. After setting up, place the plate into the instrument and run the program CNTDNA16.
- f. Prepare a mixture of 96-100% ethanol and Magnetic Beads A as per Table 3. For ease of pipetting, it is recommended to prepare a bulk mix containing 96-100% ethanol and Magnetic Beads A. We recommend a 10% overage to compensate for pipetting errors. Make sure to vortex the bulk mix before adding it to the lysate. It is recommended to prepare the bulk fresh every time a new DNA isolation is performed.
- g. After 1 hour of incubation, the instrument will pause and will allow the user to add mixture of ethanol magnetic beads prepared in step f. Open the cover and remove the plate from the instrument and add the prepared mixture from step f to Column 1 and/or Column 7.


Optional RNase A Treatment:

If RNA-free genomic DNA is required, add the equivalent of 10 KUnitz of RNase A (not to exceed 20 µL) to the lysate after addition of **96-100% ethanol** and **Magnetic Beads A** mixture to **Column 1 and/or 7** during step **f**. If RNase treatment is performed, make sure to use **CNTDNA16R** (IsoPure[™] Mini) **program instead of CNTDNA16** to include an incubation step of 37°C for 15 minutes.

- h. Put the plate back into the instrument and press OK. The instrument will continue to perform the extraction.
- i. The instrument will prompt after the procedure is complete. Remove the plate from the instrument and the DNA is now ready for further downstream processing.
- j. Transfer the eluted DNA from **Column 5 and/or Column 11** to a fresh plate. The purified DNA sample may be stored at 4°C for a few days. It is recommended that samples be placed at –20°C for long-term storage.

Flow Chart

Procedure for Norgen's Cells and Tissue DNA Isolation Kit (Magnetic Bead System) - Manual Method

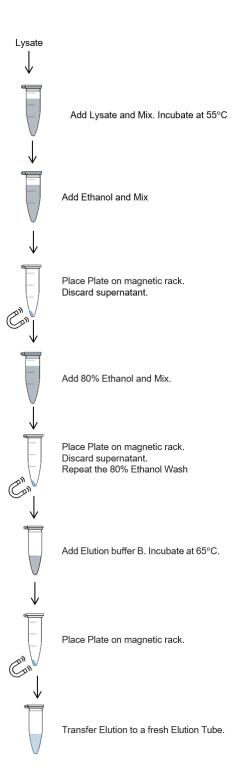

II. Manual Genomic DNA Isolation from All Types of Lysate

1A. Lysate Preparation from Animal Tissues

Notes Prior to Use

- Fresh or frozen tissues may be used for the procedure.
- Tissues should be immediately frozen and stored at -20°C or -70°C. Tissues may be stored at -70°C for several months.
- It is recommended that no more than 20 mg of tissue be used.

1A. Lysate Preparation from Animal Tissues

- a) Excise up to 20 mg of tissue sample. Either frozen or fresh tissue may be used.
- **b)** Homogenize the sample into a fine powder in liquid nitrogen using a mortar and pestle, then transfer into a nuclease-free microfuge tube.
- c) Add 300 μ L of the Lysis Buffer B to the tissue sample.

Note: The sample can also be homogenized using tools such as a rotor-stator homogenizer or a microfugesize pestle in the provided Lysis Solution.

- d) Use a syringe with a 22G needle to further homogenize the lysate by passing through the needle 10 times.
- e) Add 20 µL of Proteinase K (vortex before use) to the lysate. Mix well by gentle vortexing.
- f) Incubate at 55°C for 1 hour.

Note: Incubation times may fluctuate between 45 minutes to over 2 hours depending on the type of tissue being lysed. Lysis is considered complete when a relatively clear lysate is obtained.

g) Add 300 µL of nuclease-free water. Mix by vortexing.

Optional RNase A Treatment:

If RNA-free genomic DNA is required, add the equivalent of 10 KUnitz of RNase A (not to exceed $20 \ \mu$ L) to the lysate. Incubate at 37°C for 15 minutes.

- h) Add 600 μL of 96 100% ethanol and 40 μL of Magnetic Beads A (vortex prior to use) to the lysate collected above.
- i) Rotate the tubes on a tube revolver rotator at 10 rpm at room temperature for 5 minutes. In case the tube revolver rotator is not available, incubate the tubes at room temperature for 5 minutes and occasionally vortex the tube.
- j) Proceed to Section 2: Genomic DNA Isolation from All Types of Lysate.

1B. Lysate Preparation from Cultured Animal Cells

Notes Prior to Use

- Cells grown in suspension or monolayer may be used.
- The maximum recommended input of cells is 3 x 10⁶.
- A hemocytometer can be used in conjunction with a microscope to count the number of cells. As a general guideline, a confluent 3.5 cm plate of HeLa cells will contain 10⁶ cells.
- Cell pellets can be stored at -70°C for later use or used directly in the procedure. Determine the number of cells present before freezing.
- Frozen cell pellets should not be subjected to repeated cycles of freeze and thaw prior to beginning the protocol.

1B (i). Cell Lysate Preparation from Cells Growing in a Monolayer

- a) Detach cells by standard trypsinization method or cell scraper. Transfer an appropriate amount to a 1.5 mL microfuge tube (not provided). The maximum recommended input of cells is 3 x 10⁶.
- b) Collect cells by centrifugation at no more than 200 x g (~2,000 RPM) for 10 minutes. Discard the supernatant.
- c) Add 100 µL of phosphate buffered saline (PBS) to the cell pellet. Mix by gentle vortexing.
- d) Add 20 µL of **Proteinase K** (vortex before use) to the lysate.
- e) Add 300 µL of Lysis Buffer B. Mix well by vortexing.
- f) Incubate at 55°C for 20 minutes.

Optional RNase A Treatment:

If RNA-free genomic DNA is required, add the equivalent of 10 KUnitz of RNase A (not to exceed 20μ L) to the lysate. Mix well and incubate at 37°C for 15 minutes.

- g) Add 400 μL of **96 100% ethanol** and 40 μL of **Magnetic Beads A** (vortex prior to use) to the lysate collected above.
- h) Rotate the tubes on a tube revolver rotator at room temperature for 5 minutes. In case the tube revolver rotator is not available, incubate the tubes at 10 rpm at room temperature for 5 minutes and occasionally vortex the tube.
- i) Proceed to Section 2: Genomic DNA Isolation from All Types of Lysate.

1B (ii). Cell Lysate Preparation from Cells Growing in Suspension and Lifted Cells

- a) Transfer an appropriate amount of cells to a 1.5 mL microfuge tube (not provided). The maximum recommended input of cells is 3 x 10⁶.
- b) Collect cells by centrifugation at no more than 200 x g (~2,000 RPM) for 10 minutes. Discard the supernatant.
- c) Add 100 µL of phosphate buffered saline (PBS) to the cell pellet. Mix by gentle vortexing.
- d) Add 20 μ L of **Proteinase K** (vortex before use) to the lysate.
- e) Add 300 µL of Lysis Buffer B. Mix well by vortexing.
- f) Incubate at 55°C for 20 minutes.

Optional RNase A Treatment:

If RNA-free genomic DNA is required, add the equivalent of 10 KUnitz of RNase A (not to exceed 20 μ L) to the lysate. Mix well and incubate at 37°C for 15 minutes.

- g) Add 400 μL of **96 100% ethanol** and 40 μL of **Magnetic Beads A** (vortex prior to use) to the lysate collected above.
- h) Rotate the tubes on a tube revolver rotator at 10 rpm at room temperature for 5 minutes. In case the tube revolver rotator is not available, incubate the tubes at room temperature for 5 minutes and occasionally vortex the tube.
- i) Proceed to Section 2: Genomic DNA Isolation from All Types of Lysate.

1C. Lysate Preparation for Viral DNA

Notes Prior to Use

- For the isolation of integrated viral DNA, follow **Section 1A** if the starting material is animal tissue and follow **Section 1B** if the starting material is cell culture.
- For the isolation of DNA from free viral particles, follow the protocol provided below.
- Up to 150 µL of viral suspension can be processed.
- Fresh samples are recommended. Frozen samples may be used, however the yield of genomic DNA may be decreased.

1C. Lysate Preparation for Viral DNA

- a) To a 1.5 mL microfuge tube (not provided), add 20 µL of Proteinase K (vortex before use).
- b) Transfer up to 150 μL of viral suspension to the tube.
- c) Add 300 µL of Lysis Buffer B. Mix well by vortexing.
- d) Incubate at 55°C for 1 hour.

Note: Incubation times may fluctuate between 45 minutes to over 2 hours depending on the type of cell being lysed. Lysis is considered complete when a relatively clear lysate is obtained.

Optional RNase A Treatment:

If RNA-free genomic DNA is required, add the equivalent of 10 KUnitz of RNase A (not to exceed $20 \ \mu$ L) to the lysate. Mix well and incubate at 37°C for 15 minutes.

- e) Add 450 μL of 96 100% ethanol and 40 μL of Magnetic Beads A (vortex prior to use) to the lysate collected above.
- f) Rotate the tubes on a tube revolver rotator at 10 rpm at room temperature for 5 minutes. In case the tube revolver rotator is not available, incubate the tubes at room temperature for 5 minutes and occasionally vortex the tube.
- g) Proceed to Section 2: Genomic DNA Isolation from All Types of Lysate.

Section 2. Genomic DNA Isolation from All Types of Lysate

Note: The remaining steps of the procedure for the purification of genomic DNA are the same from this point forward for all the different types of lysate.

- a) Assemble a magnetic separation rack and place the sample tube in the magnetic rack. Allow to sit for 1 minute. Incubate for more time if the solution is not clear.
- b) Aspirate and discard supernatant without touching the magnetic beads.
- c) Remove the sample tube from the magnetic rack and gently add 500 μL of **Solution WN (ensure ethanol was added).** Resuspend by vortexing or pipetting and incubate at room temperature for 1 minute.
- d) Place the sample tube on the magnetic rack and allow to sit for 1 minute. Incubate for more time if the solution is not clear.
- e) Aspirate and discard supernatant without touching the magnetic beads.
- f) Remove the sample tube from the magnetic rack and gently add 500 μL of freshly prepared **70% ethanol**. Resuspend by vortexing or pipetting and incubate at room temperature for 1 minute.
- g) Place the sample tube on the magnetic rack and allow to sit for 1 minute. Incubate for more time if the solution is not clear.
- h) Aspirate and discard supernatant without touching the magnetic beads.
- i) Repeat **Steps 2f 2h** for a second wash step.

Note: Remove as much of the 70% ethanol in the sample tube as possible by pipetting.

- j) Incubate the open tube at 65°C for 5 minutes to dry the magnetic beads.
- k) Remove the sample tube from the magnetic rack and add 50-100 μL of **Elution Buffer B**. Mix by vortexing and incubate at 65°C for 10 minutes.
- I) Briefly vortex and place sample tube on the magnetic rack and allow to sit for 1 minute. Incubate for more time if the solution is not clear.
- m) Carefully transfer the elution to a fresh 1.7 mL elution tube (provided) without touching the magnetic beads.
- n) The purified DNA sample may be stored at 4°C for a few days. It is recommended that samples be placed at -20°C for long-term storage.

Troubleshooting Guide

Problem	Possible Cause	Solution and Explanation	
Magnetic beads were accidentally pipetted up with the supernatant.	The pipette tip was placed too close to the magnetic beads while pipetting	Return the magnetic beads and the supernatant back into the sample tube. Mix well, and place the tube back onto the magnetic separation rack for the specified time. Carefully remove the supernatant without touching the magnetic beads.	
	Incomplete lysis of cells	Ensure that correct lysis protocol was applied to the sample. Ensure Proteinase K is added properly.	
The yield of genomic DNA is	Amount of magnetic beads added was not sufficient	Ensure that the Magnetic Bead A is mixed well prior to use to avoid any inconsistency in DNA isolation.	
low	DNA concentration in the cell or tissue sample being used is low.	Some samples contain very little target DNA. This varies from individual to individual based on numerous variables. Extend the incubation time of Proteinase K digestion or reduce the amount of tissue or cells used for lysis.	
Very gelatinous prior to adding	The lysate solution mixture is not homogeneous	To ensure a homogeneous solution, vortex for 10-15 seconds before adding the magnetic beads to the lysate.	
the Magnetic bead and Ethanol	Maximum number of cells or amount of tissue exceeds kit specifications	Refer to specifications to determine if amount of starting material falls within kit specifications.	
DNA does not perform well in downstream	DNA was not washed with 70% Ethanol	Traces of salt from the binding step may remain in the sample if the magnetic beads are not washed with 70% Ethanol. Salt may interfere with downstream applications, and thus must be washed from the magnetic beads.	
applications.	Ethanol carryover	Ensure that the drying step after the 70% ethanol wash steps is performed, in order to remove traces of ethanol prior to elution. Ethanol is known to interfere with many downstream applications.	
RNA is present in eluted DNA.	RNA is co-eluted with the DNA	Carry out a digestion with RNase A on the elution if the RNAse present will interfere with downstream applications. Refer to manufacturer's instructions regarding amount of enzyme to use, optimal incubation time and temperature.	

Related Products	Product #
Cells and Tissue DNA Isolation 96-Well Kit (Magnetic Bead System)	62500
Cells and Tissue DNA Isolation Kit (50 Prep)	53100
Bacterial Genomic DNA Isolation Kit	17900
Blood DNA Isolation Mini Kit	46300
Adhesive Tape – 100 Units	28394

Technical Support

Contact our Technical Support Team between the hours of 9:00 and 5:30 (Eastern Standard Time) at (905) 227-8848 or Toll Free at 1-866-667-4362. Technical support can also be obtained from our website (www.norgenbiotek.com) or through email at techsupport@norgenbiotek.com.

3430 Schmon Parkway, Thorold, ON Canada L2V 4Y6 Phone: (905) 227-8848 Fax: (905) 227-1061 Toll Free in North America: 1-866-667-4362

PI59100-4